Powertube®

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Seite 3

Mode d'emploi

page 15

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pagina 29

Instruction Manual

page 41

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
Instruction manual

This instruction manual is considered an integral part of the device. It shall be kept within reach at all times. Exact compliance with these instructions is essential for the correct handling of the device as well as for the safety of patient and user.


Tabel of Contents Page

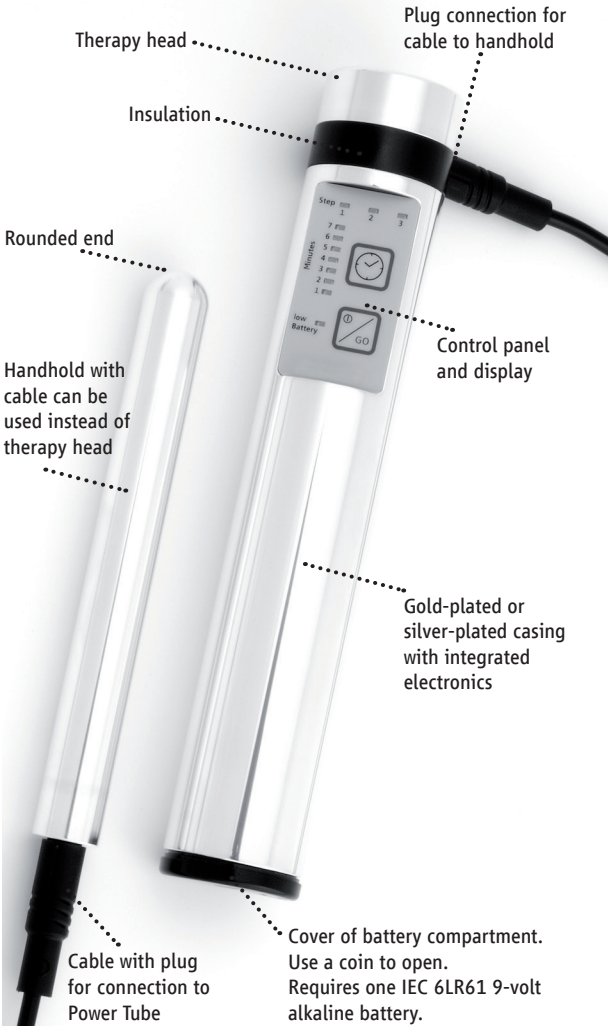1.	Illustration of the device	42
2.	Control panel with functions	43
3.	Description of the device	44
4.	Quick start guide	44
5.	Battery and battery compartment	45
6.	Specifications	46
7.	Handling	47
8.	Indications	48
9.	Contraindications	48
10.	Safety tips	48
11.	Care, maintenance and storage	49
12.	Errors and malfunctions, function test	50
13.	Accessories	50
14.	Guarantee	50
15.	Safety inspection	51
16.	CE-Conformity declaration	52

1. Illustration of the device

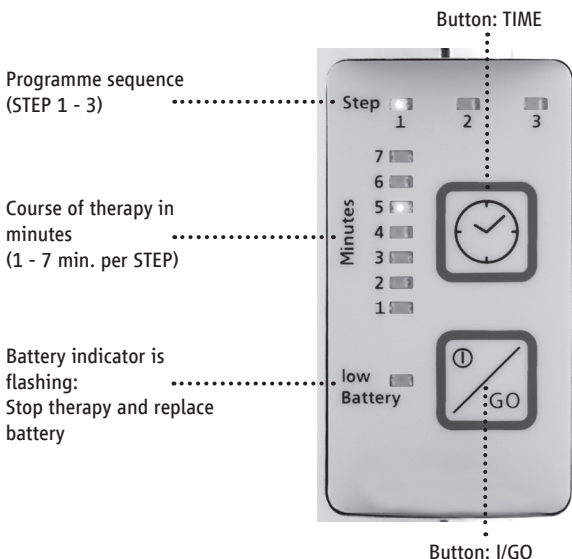2. Control panel with functions

-  Turning on the device (push button briefly)
-  Select number of minutes (1-7 minutes)
-  Start therapy (push button briefly) (manual STEP selection)
-  Turning off the device (push button for an extended period of time or device will automatically turn off after 30 seconds)

3. Description of the device

Intended use

The POWER TUBE is intended for the electronic stimulation of nerves. It is used for self-treatment of symptoms and pain.

Its application does not replace treatment or instructions from a doctor, but has a positive influence on symptoms and pain. Using the device can alleviate pain.

Function designation

The transcutaneous electronic nerve stimulation device POWER TUBE is an easy to use, battery-operated electronic device designed to stimulate the nerves.

Corresponding nerve points are stimulated directly with electronic impulses. The device can be operated by persons who do not have formal medical training.

The device is described in the enclosed easy to follow instruction manual to make self-treatment a straightforward process. The handling of the device is intended to suit users with no medical training. Various applications as well as corresponding therapy times are outlined in the instruction manual.

4. Quick start guide

Before using for the first time

Immediately after receipt, the device should be inspected for transport damage and be subject to a function test. In case of transport damage, the shipper shall be notified immediately and a damage report be drawn up to support any claim for damages. In addition to name and address of the recipient, the protocol shall include the date of receipt, type and serial number of the device as well as a description of the damage. The original packaging should be kept during the guarantee period to ensure that the device can be returned in its original packaging if necessary.

Every device is tested by the manufacturer in terms of func-


tionality and safety before being shipped. Please observe the following to ensure that the device continues to operate properly following transport and installation:

The operator should only start using the device after carrying out a function test (see chapter 12).

5. Battery and battery compartment

Use alkaline batteries only and check the polarities (plus/minus) carefully. By using this kind of battery you may be certain to end your therapy time even in case the battery is almost empty by showing the sign «low battery». The POWER-TUBE device is adjusted to this type of batteries.

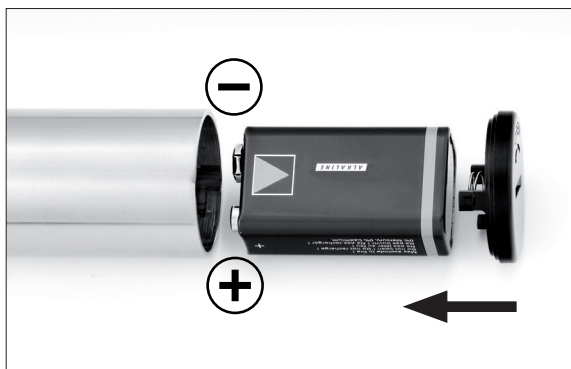If you use rechargeable accus (e.g. NiMH), the proper functioning of the sign «low battery» is not completely reliable because of the different mode of discharging. The POWER-TUBE device may abruptly stop functioning without prior warning. If you still prefer to use rechargeable accus, please recharge them periodically.

6. Specifications

Display:	Time of therapy in minutes, 7 light-emitting diodes Therapy steps 1 - 3 Exchange light-emitting diode battery
Error messages:	Intermittent tone: handholds or contact areas are not applied properly Continuous tone: handholds have short circuited
Battery:	IEC 6LR61 9-volt block battery
Output voltage:	AC, maximum 9-volt pp and minimum 6-volt pp
Duration of treatment:	3 - 21 minutes
Programme sequence:	3 steps
Measurements:	Length: 18 cm, diameter: 3.1 cm
Weight with softbag:	≈ 750 g
Application class:	BF 
Safety class:	Ila unit
Operation:	10°C - 40°C Humidity 30% - 75% without condensation, air pressure 700 hPa - 1060 hPa
Please refer to instruction manual!	
Production Date	month / year-serial number

7. Handling

1. Turning on the device: push I/GO button briefly.
2. Set appropriate therapy time per step by means of TIME button. Start with one minute per step, i.e. the complete treatment time is 3 minutes. Increase the time with every treatment by 1 minute per step. Complete steps 1 and 2 (see below and beside) once a day.
3. Place device on trigger point.
4. Start treatment by pushing the I/GO button.
5. You will hear a beep signal after completion of every step.
6. The end of the treatment is announced with two beep signals.
7. Turning off the device: push I/GO button for a prolonged period of time or device turns off automatically after 30 seconds.

Error message: Insufficient or no skin contact is indicated with continuous beep signals.


Step 1

Local treatment at trigger point

Step 2

Treatment of abdominal area

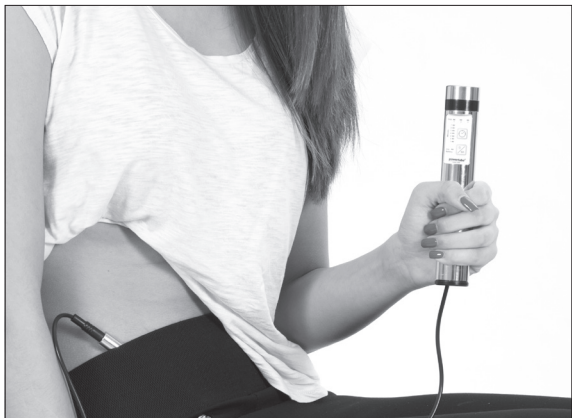Please ensure good body contact. Moisten skin with water if necessary.

Consume large quantities of water before and after the treatment.

8. Indications

- Pain relief

9. Contraindications

- Pregnancy: do not use the device.
- Patients with a cardiac pacemaker: do not use the device.
- Open wounds: do not use the device.
- Fresh scars: do not use the device.

10. Safety tips

- POWERTUBE may only be used for the electronic stimulation of nerves and musculature.
- Do not use the device for patients who have been fitted with a cardiac pacemaker.

- Pay attention to strangulation by cables, especially with excess lengths.
- Therapy time for the first application (without cables, direct therapy)
 - Head max. 2-2-2 minutes
 - Chest max. 3-3-3 minutes
- After a week, the therapy time can be increased slightly.
- Simultaneous connection with other devices is not permitted.
- POWERTUBE may only be operated with a 9-volt battery.
- POWERTUBE may only be operated with the accessories included in the range of delivery.
- POWERTUBE should be disposed of as hazardous waste.
- For technical questions, handling or maintenance POWERTUBE, please contact the manufacturer.

Swisspowertube AG is only responsible for the reliability and safety of POWERTUBE if:

» **POWERTUBE is operated in accordance with the instruction manual**

» **expansions and modifications are carried out only by persons who are authorised by Swisspowertube AG**

11. Care, maintenance and storage

Effective protection of the POWERTUBE from damage includes:

- Appropriate operation and maintenance.
- Protection against moisture, high temperatures, dust and direct sunlight.
- Protection from children, animals and pests.
- After each use, clean the contact surfaces of the POWERTUBE device and the hand electrodes with a commercially available cleaning agent (e.g. 1-propanil, ethanol, glyoxal). The print nameplate should be left out.
- If the device is not used for a longer period of time, the remove the battery from the device.
- Dispose of the battery as hazardous waste.
- The expected service life of new alkaline batteries is approx. 2 months (with daily use of 20 minutes).

12. Errors and malfunctions, function test

- When the battery is discharged, this is indicated by a battery indicator LED.
- After switching on, POWERTUBE automatically performs a battery test.
- If the hand electrode or the POWERTUBE housing is not correctly applied during therapy, a warning signal will sound.

Test the proper functioning of your POWERTUBE device yourself

1. Insert a new or newly charged battery into the battery compartment.
2. Connect the handhold with the cable for connection to the POWERTUBE device.
3. Push the start-button on the control panel.
4. Let the standard setting of 1 minute therapy (basic setting).
5. Start the therapy by pushing the start-button again and touch with the handhold the POWERTUBE casing (not the therapy head!).
6. If the POWERTUBE device now emits a continuing peep-sound, the device functions perfectly.
7. Switch off the POWERTUBE and start over with the therapy.

13. Accessories

- Hand electrodes
- Instruction manual
- IEC 6LR61 9-volt battery
- Portable case for POWERTUBE

14. Guarantee

- Device and accessories shall be inspected for transport damage immediately after receipt. Corresponding claims for damages can only be asserted if the vendor or shipper have been notified immediately. A damage report shall be compiled.
- The guarantee for POWERTUBE is two years from the day of delivery.
- The extent of guarantee includes free repair of the device, provided the damage was caused by a material or manufacturing

error.

- Additional claims, particularly claims for damages, are excluded.
- The repair may only be carried out by Swisspowertube AG, one of its agencies or authorised specialist dealers.
- The guarantee shall become void if improper modifications or repairs are carried out.
- The guarantee period will neither be extended nor renewed because of guarantee work.
- Swisspowertube AG reserves all rights to this instruction manual, particularly the right of duplication and distribution as well as translations. No part of this instruction manual may be copied or distributed in any form without prior written consent from Swisspowertube AG.
- Informations contained in this instruction manual may be modified or upgraded by Swisspowertube AG.
- Updates of Manuals are provided for downloading on our Website www.swisspowertube.ch.

15. Safety inspection

- § 6 of the directive on installing, operating and using active medical products stipulates regular safety inspections to guarantee the safety of the device against age, wear and tear etc.
- The operator is responsible for having these safety inspections carried out on the device to the prescribed extent and in due time. Safety inspections may only be carried out by the manufacturer or by explicitly authorised persons.


The following safety inspections are required for POWERTUBE:


- Examination of labelling and instruction manual
- Visual inspection of device and accessories for damages
- Function test of all sensors and control lamps on the device

16. EG-Conformity declaration

CE 0482

Name and address of the company	Swisspowertube AG Breite 9 CH-3636 Forst Switzerland SRN-Nr: CHRN-MF-20000238
We declare in the position as solely responsible company that	
the medical device	POWERTUBE (gold) Typ FPT.B1 POWERTUBE (silver) Typ FPT-D.B2
of class	Ila rule 9 according to Annex IX of Directive 93/42/EEC
meets the provisions of the directive 93/42/EEC and its transpositions in national laws which apply to it. The declaration relates to the medical products supplied.	
Conformity assessment procedure	Directive 93/42/EEC Annex II without section 4
Notified body	MEDCERT GmbH Pilatuspool 2 D-20355 Hamburg Deutschland  0482
CH-Forst, 5th of August 2021	Martin Frischknecht, CEO 
Name and address of the EU authorized representative	Wellbalance GmbH Klosterstrasse 10 D-77716 Haslach Germany SRN-Nr: DE-AR-000010274

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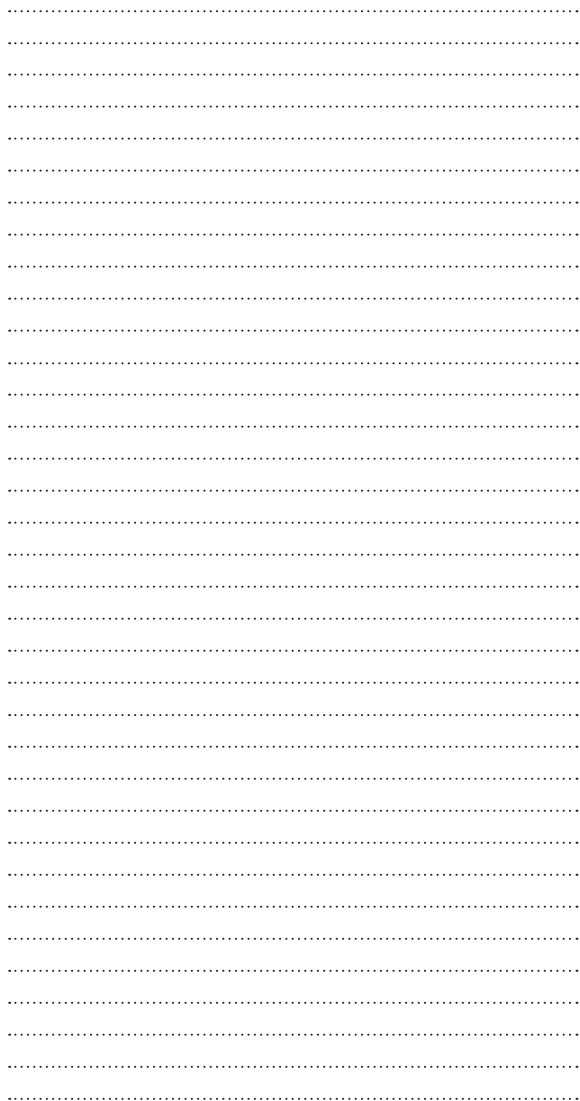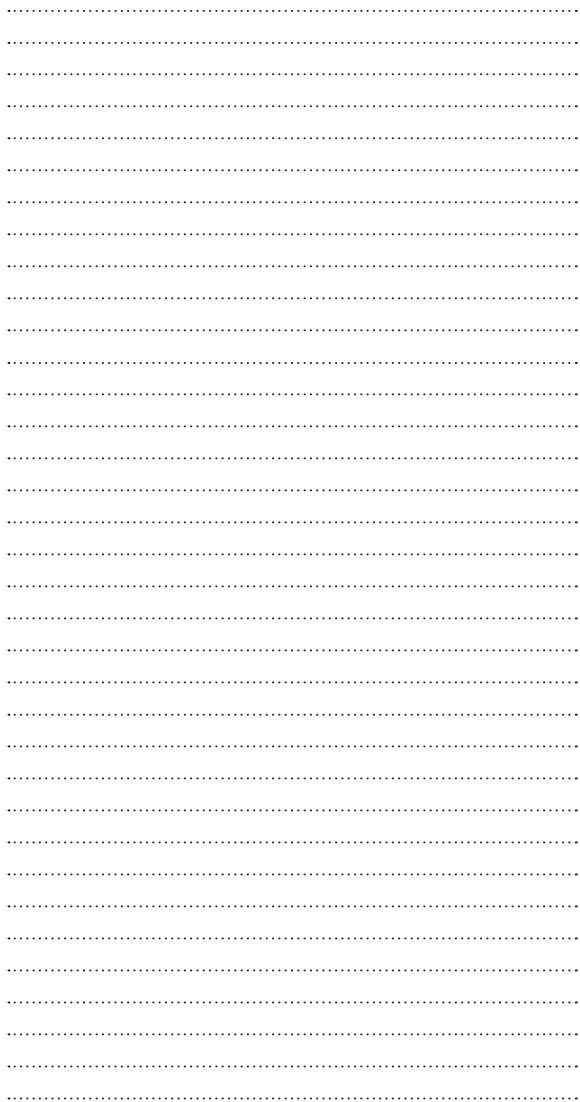©2021 Swisspowertube AG. All rights reserved. All trade names are registered trademarks of

Swisspowertube AG
Breite 9
3636 Forst
Switzerland

Edition 2021-08-05
Printed in Switzerland

Powertube®